



Applied Neuroscience, Inc.

3831 Tyrone Blvd. N., Suite 103, Largo, FL 33709 USA
 Phone: (727) 657-4349 Email: qeeg@appliedneuroscience.com

Quantitative EEG Analyses

PATIENT INFORMATION

Name: John Doe
 EEG ID: EEG01
 Age: 54.61
 Gender: Male
 Handedness: Right
 Eyes: Closed

RECORDING

Date: 12/12/2003
 Ref. By:
 Test Site:
 Analysis Length: 02:19
 Ave. SH Reliability: 0.99
 Ave. TRT Reliability: 0.98

MEDICATION: Zoloft, Dilantin

HISTORY: He was hit on the right side of the head with a baseball bat. No immediate LOC. Next day he was found prone on the floor, delirious, with difficulty responding to verbal commands. He was taken to the ER and CT Scan revealed a right hemisphere epidural hematoma with pocket of blood in the right occipital region.

QEEG Impression

SUMMARY: The qEEG analyses were abnormal and showed dysregulation in bilateral frontal lobes - especially in the right frontal lobe, bilateral temporal lobes - especially in the right temporal lobe, bilateral parietal lobes - especially in the right parietal lobe, and bilateral occipital lobes - especially in the right occipital lobe. LORETA showed dysregulation in the right supramarginal gyrus, right postcentral gyrus, and right posterior cingulate & cuneus. The frontal lobes are involved in executive functioning, abstract thinking, expressive language, sequential planning, mood control, and social skills. The temporal lobes are involved in auditory information processing, short-term memory, receptive language on the left, and face recognition on the right. The parietal lobes are involved in visual-spatial information processing, short-term memory, executive attention, receptive language on the left, and empathy control and awareness of emotional expression in others on the right (e.g., prosody). The occipital lobes are involved in the visual processing of color, form, movement, visual perception, and spatial processing. The postcentral gyrus is involved with the sense of touch and efferent motor feedback. The posterior cingulate gyrus is involved in spatial memory, autobiographical memory, salience and attention control by regulation of limbic emotional and memory input to the hippocampus and neocortex. The cuneus is involved in basic visual processing and is modulated by extraretinal effects, like attention, working memory, and reward expectation. To the extent there is deviation from normal electrical patterns in these structures, then sub-optimal functioning is expected.

Robert W. Thatcher

Dr. Robert W. Thatcher, PhD, QEEG-D

Detailed Narrative

LINKED EARS: The Linked Ears power spectral analyses were deviant from normal with excessive power in bilateral frontal regions - especially in the right frontal region over a wide frequency range. Excessive power was present in bilateral temporal regions - especially in the right temporal region over a wide frequency range. Excessive power was present in bilateral parietal regions - especially in the right parietal region over a wide frequency range. Excessive power was also present in bilateral occipital regions - especially in the right occipital region from 1 - 2 Hz.

SURFACE LAPLACIAN: The Laplacian power spectral analyses were deviant from normal with excessive power in bilateral frontal regions - especially in the right frontal region over a wide frequency range. Excessive power was present in bilateral temporal regions - especially in the right temporal region at 1 Hz. Excessive power was present in the right parietal region over a wide frequency range. Excessive power was also present in bilateral occipital regions - especially in the right occipital region over a wide frequency range.

NEUROIMAGING: LORETA 3-dimensional source analyses were consistent with the surface EEG and showed elevated current sources in the right posterior cingulate & cuneus with a maximum at 1 Hz (Brodmann Area 30). Elevated LORETA current source also were present in the right supramarginal gyrus with a maximum at 2 Hz (Brodmann Area 7). Elevated LORETA current source was present in the right supramarginal gyrus with a maximum at 4 Hz (Brodmann Area 7). Elevated LORETA current source was present in the right supramarginal gyrus with a maximum at 5 Hz (Brodmann Area 7). Elevated LORETA current source was present in the right postcentral gyrus with a maximum at 6 Hz (Brodmann Area 1).

CONNECTIVITY ANALYSES: EEG amplitude asymmetry, coherence and EEG phase were deviant from normal, especially in frontal, temporal, parietal and occipital relations. Elevated coherence was present in frontal, temporal, parietal and occipital regions which indicates reduced functional differentiation. Reduced coherence was present in frontal, temporal, parietal and occipital regions which indicates reduced functional connectivity. Both conditions are often related to reduced speed and efficiency of information processing.

DISCRIMINANT ANALYSES: The mild head injury discriminant function detected a pattern in the EEG that is commonly present in individuals with a history of mild traumatic brain injury.

General Neurofeedback Recommendations

Surface Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta		Theta		Alpha		Beta1		Beta2		Beta3			High Beta																	
	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30	
FP1																															
FP2																															
F3																															
F4																															
C3																															
C4																															
P3																															
P4																															
O1																															
O2																															
F7																															
F8																															
T3																															
T4																															
T5																															
T6																															
Fz																															
Cz																															
Pz																															

Coherence Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta		Theta		Alpha		Beta1		Beta2		Beta3			High Beta																	
	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30	
F8-T3																															
F8-T5																															
C4-Cz																															
O1-F8																															
P3-F8																															
P4-Cz																															
C4-Fz																															
C3-P4																															
T3-T4																															
C3-C4																															

Loreta Neurofeedback: (move towards Z = 0)

FREQ BAND	L	Delta		Theta		Alpha		Beta1		Beta2		Beta3			High Beta																
		1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
BA7	R																														
BA5	R																														
BA1	R																														
BA4	R																														
BA40	R																														
BA3	R																														
BA2	R																														
BA30	R																														
BA35	R																														
BA34	R																														
BA36	R																														
BA31	R																														
BA28	R																														
BA27	R																														
BA19	R																														
BA37	R																														
BA41	R																														
BA20	R																														
BA18	R																														
BA29	R																														

The implications for neurotherapy are offered based upon the clinical evaluation of the patient as well as the reference data base results. These suggestions for neurotherapy should be evaluated with caution and should only be considered as possible strategies that the clinician may have considered in his/her evaluation. If the patient is depressed, then the clinician should consider treating this condition first through alpha frequency enhancement or some other biofeedback protocol that may reduce depression. If depression or poor mood that may reduce depression. If depression or poor mood that may reduce depression. If depression or poor mood that may reduce depression. If depression or poor mood that may reduce depression. If depression or poor mood that may reduce depression. If depression or poor mood that may reduce depression. If depression or poor mood that may reduce depression. If depression or poor mood that may reduce depression.

Addiction/Reward Network Recommendations

Surface Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta				Alpha				Beta1			Beta2			Beta3			High Beta										
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30	
F3	■																														
F4	■																														
C3	■																														
C4	■																														
F7																															
F8	■																														
Fz	■																														
Cz	■																														

Coherence Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta				Alpha				Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30		
C4-Cz				■																												
C4-Fz				■																												
C3-C4				■																												
F4-C4				■																■												
C3-F8	■																															
F3-F8	■																															
F3-C4				■																												
F8-Cz	■																															
C3-F7				■																												
F8-Fz	■																															

Loreta Neurofeedback: (move towards Z = 0)

FREQ BAND	L	Delta			Theta				Alpha				Beta1			Beta2			Beta3			High Beta										
FREQ (Hz)	R	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30	
BA34	R	■																														
BA34	L	■																														
BA13a	L	■																														
BA47	L	■																														
BA25	R	■																														
BA25	L	■																														
BA13a	R	■																														
BA45	L	■																														
BA47	R	■																														
BA46	L	■																														
BA44	R	■																														
BA44	L	■																														
BA24	R	■																														
BA24	L	■																														
BA45	R	■																														
BA46	R	■				■																										
BA32	R	■																														
BA32	L	■																														

Anxiety Network Recommendations

Surface Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta																			
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30								
F3																																						
F4																																						
C3																																						
C4																																						
Fz																																						
Cz																																						

Coherence Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta																						
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30											
C4-Cz																																									
C4-Fz																																									
C3-C4																																									
F4-C4																																									
F3-C4																																									

Loreta Neurofeedback: (move towards Z = 0)

FREQ BAND	L	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta																			
FREQ (Hz)	R	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30								
BA7	R																																						
BA4	R																																						
BA13a	L																																						
BA21	L																																						
BA6	R																																						
BA13p	L																																						
BA13p	R																																						
BA13a	R																																						
BA10	L																																						
BA21	R																																						
BA10	R																																						

Attention Network Recommendations

Surface Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
FP1																														
FP2																														
F3		█																												
F4		█	█																											
P3		█																												
P4		█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
Fz		█																												
Cz		█																												
Pz		█																												

Coherence Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
P4-Cz	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
P3-P4				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
P4-Fz				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
F4-P4				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
FP1-P3				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
P4-Pz				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
F3-P3	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
FP2-P3				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
P3-Fz	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
F3-P4				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█

Loreta Neurofeedback: (move towards Z = 0)

FREQ BAND	L	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	R	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
BA7	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA40	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA19	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA39	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA6	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA8	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA9	L	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA19	L	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	

Default Mode Network Recommendations

Surface Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
F3	█																													
F4	█	█																												
P3	█																													
P4	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
O1	█																													
O2	█	█																												
T5	█																													
T6	█	█																												
Fz	█																													
Cz	█																													
Pz	█																													

Coherence Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
P4-Cz	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
P3-P4				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
F3-T5	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
O2-T6																														
P4-Fz				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
F3-O1	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
F4-P4				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
P4-Pz				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
F3-P3	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
P3-Fz	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█

Loreta Neurofeedback: (move towards Z = 0)

FREQ BAND	L	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	R	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
BA7	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA40	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA2	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA30	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA35	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA31	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA19	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA29	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA39	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA35	L	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA30	L	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA11	L	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA11	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA31	L	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA10	L	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA10	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA29	L	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA19	L	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	

Executive Function Network Recommendations

Surface Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
F3	█																													
F4	█	█																												
C3	█																													
C4	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
P3	█																													
P4	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
F7																														
F8	█																													
Fz	█																													
Cz	█																													
Pz	█																													

Coherence Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
C4-Cz				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
P3-F8	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
P4-Cz	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
C4-Fz				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
C3-P4				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
C3-C4				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
F4-C4				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
C3-F8	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
P3-P4				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█
P4-Fz				█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█

Loreta Neurofeedback: (move towards Z = 0)

FREQ BAND	L	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	R	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
BA7	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA40	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA19	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA37	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA11	L	█																													
BA11	R	█																													
BA22	L	█																													
BA22	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA37	L	█																													
BA10	L	█																													
BA10	R	█																													
BA46	L	█																													
BA8	R	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	█	
BA9	L	█																													
BA19	L	█																													
BA46	R	█	█																												

Memory Network Recommendations

Surface Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta																	
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30						
T3																																				
T4																																				
T5																																				
T6																																				

Coherence Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta																	
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30						
T3-T4																																				
T3-T6																																				
T4-T5																																				

Loreta Neurofeedback: (move towards Z = 0)

FREQ BAND	L	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta														
FREQ (Hz)	R	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30			
BA35	R																																	
BA34	R																																	
BA36	R																																	
BA28	R																																	
BA35	L																																	
BA34	L																																	
BA28	L																																	
BA36	L																																	

Pain Network Recommendations

Surface Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta																			
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30								
F4	█																																					
C4	█																																					
P3	█																																					
P4	█																																					
Fz	█																																					
Cz	█																																					
Pz	█																																					

Coherence Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta																																																											
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30																																																
C4-Cz																															█																																															
P4-Cz	█																																																																													
C4-Fz																															█																																															
F4-C4																															█																																															
P3-P4																															█																																															
P4-Fz																															█																																															
F4-P4																															█																																															
C4-Pz																															█																																															
P4-Pz																															█			█																																												
C4-P4																															█																																															

Loreta Neurofeedback: (move towards Z = 0)

FREQ BAND	L	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta																			
FREQ (Hz)	R	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30								
BA5	R	█																																					
BA1	R	█																																					
BA4	R	█																																					
BA3	R	█																																					
BA2	R	█																																					
BA13a	L	█																																					
BA13a	R	█																																					
BA5	L	█																																					
BA33	R	█																																					
BA33	L	█																																					
BA24	R	█																																					
BA24	L	█																																					
BA32	R	█																																					
BA32	L	█																																					

Mirror Neuron Network Recommendations

Surface Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta				Alpha					Beta1			Beta2			Beta3			High Beta																
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30								
C3	█																																					
C4	█																																					
F7																																						
F8	█																																					
T3	█																																					
T4	█																																					
Fz	█																																					
Cz	█																																					

Coherence Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta				Alpha					Beta1			Beta2			Beta3			High Beta																																														
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30																																						
F8-T3	█												█																																																							
C4-Cz																															█																																					
C4-Fz																															█																																					
T3-T4	█																																																																			
C3-C4																															█																																					
C3-F8	█																																																																			
F8-Cz	█																																																																			
C3-F7																															█																																					
T3-Fz	█																																	█																																		
F8-Fz	█																																																																			

Loreta Neurofeedback: (move towards Z = 0)

FREQ BAND	L	Delta			Theta				Alpha					Beta1			Beta2			Beta3			High Beta																
FREQ (Hz)	R	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30								
BA5	R	█																																					
BA1	R	█																																					
BA40	R	█																																					
BA27	R	█																																					
BA13a	L	█																																					
BA6	R	█																																					
BA27	L	█																																					
BA13a	R	█																																					
BA45	L	█																																					
BA46	L	█																																					
BA44	R	█																																					
BA5	L	█																																					
BA44	L	█																																					
BA45	R	█																																					
BA46	R	█		█																																			

Salience Network Recommendations

Surface Neurofeedback: (move towards Z = 0)

FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
FP1																														
FP2																														
F3		■																												
F4		■	■																											
P3		■																												
P4		■	■	■	■	■	■	■	■	■	■																			
Fz		■																												
Pz		■																												

Coherence Neurofeedback: (move towards Z = 0)











FREQ BAND	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
P3-P4				■	■	■	■	■	■	■	■																			
P4-Fz				■	■	■	■	■	■	■	■																			
F4-P4				■	■	■	■	■	■	■	■																			
FP1-P3				■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■
P4-Pz				■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■
F3-P3	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■
FP2-P3	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■
P3-Fz	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■
F3-P4				■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■

Loreta Neurofeedback: (move towards Z = 0)

FREQ BAND	L	Delta			Theta			Alpha			Beta1			Beta2			Beta3			High Beta											
FREQ (Hz)	R	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30
BA30	R	■	■	■																											
BA31	R	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	
BA29	R	■	■	■				■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	
BA30	L	■	■	■																											
BA23	R	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	
BA13a	L	■	■	■																											
BA25	R	■	■	■																											
BA25	L	■	■	■																											
BA22	L	■	■	■																											
BA31	L	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	
BA22	R	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	
BA13a	R	■	■	■																											
BA10	L	■	■	■																											
BA23	L	■	■	■																											
BA10	R	■	■	■																											
BA29	L	■	■	■																											
BA33	R	■	■	■																											
BA8	R	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	■	
BA33	L	■	■	■																											
BA24	R	■	■	■																											


Significant Functioning Network Z-Scores

(Circles = zSCOREs, Radial Lines = Brodmann Areas)

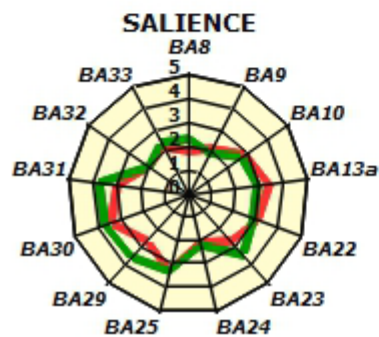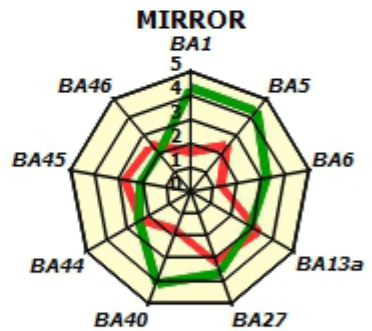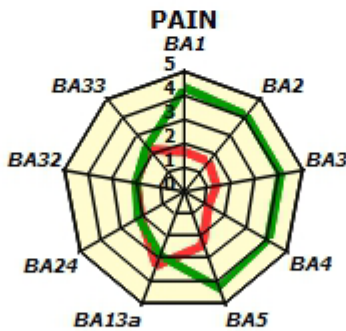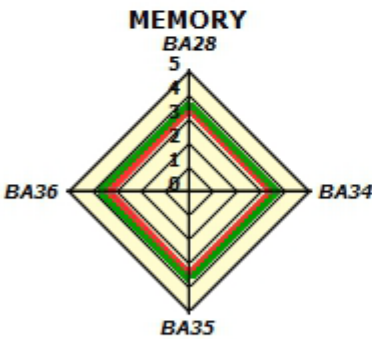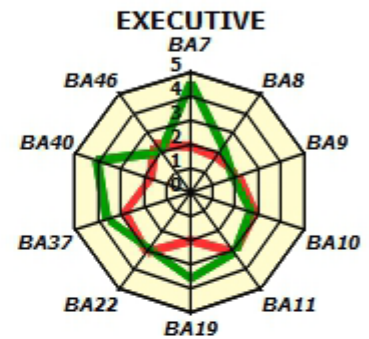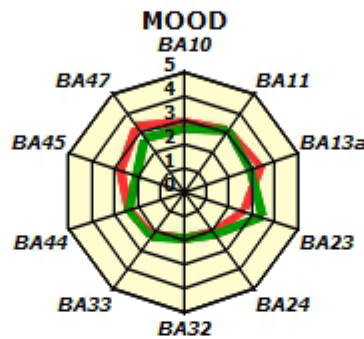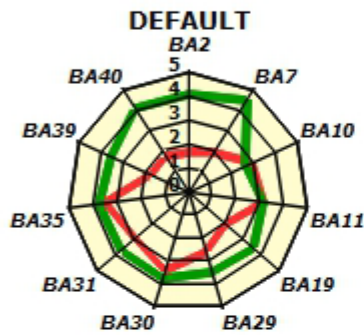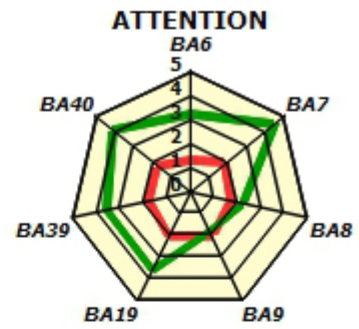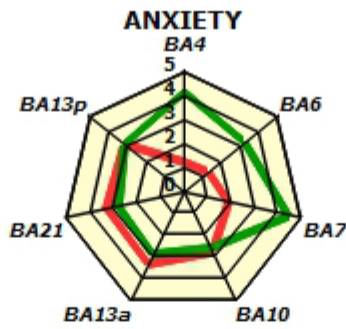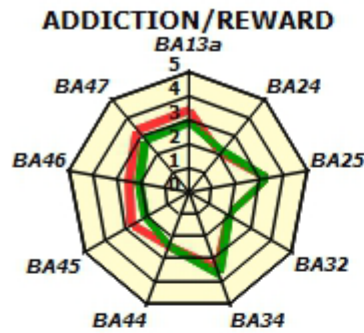


■ LEFT ■ RIGHT

QEEG Overview Maps

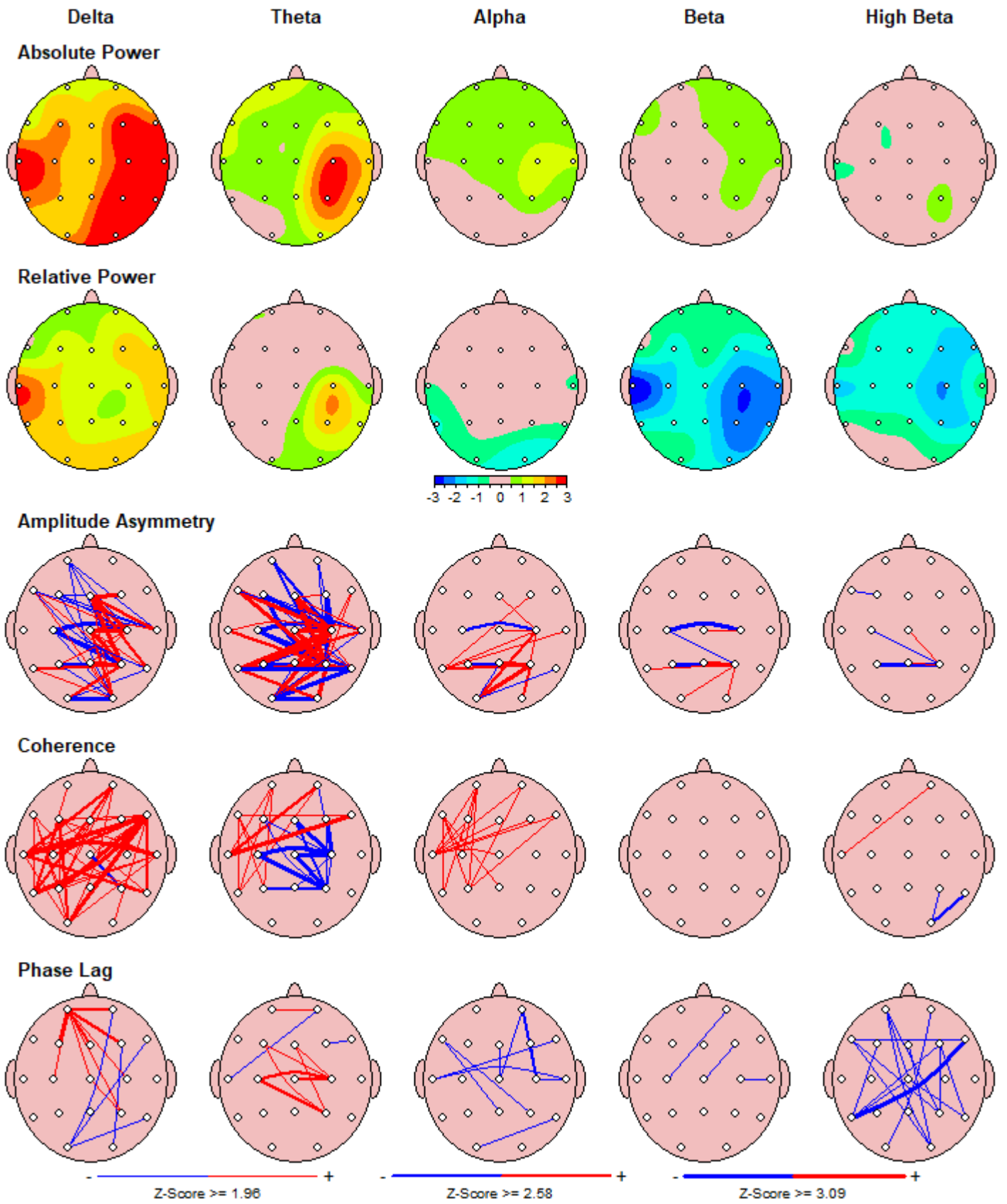


Conventional EEG Samples and Quantitative EEG Analyses

Example of Linked Ears EEG and Absolute Power - Eyes Closed Condition

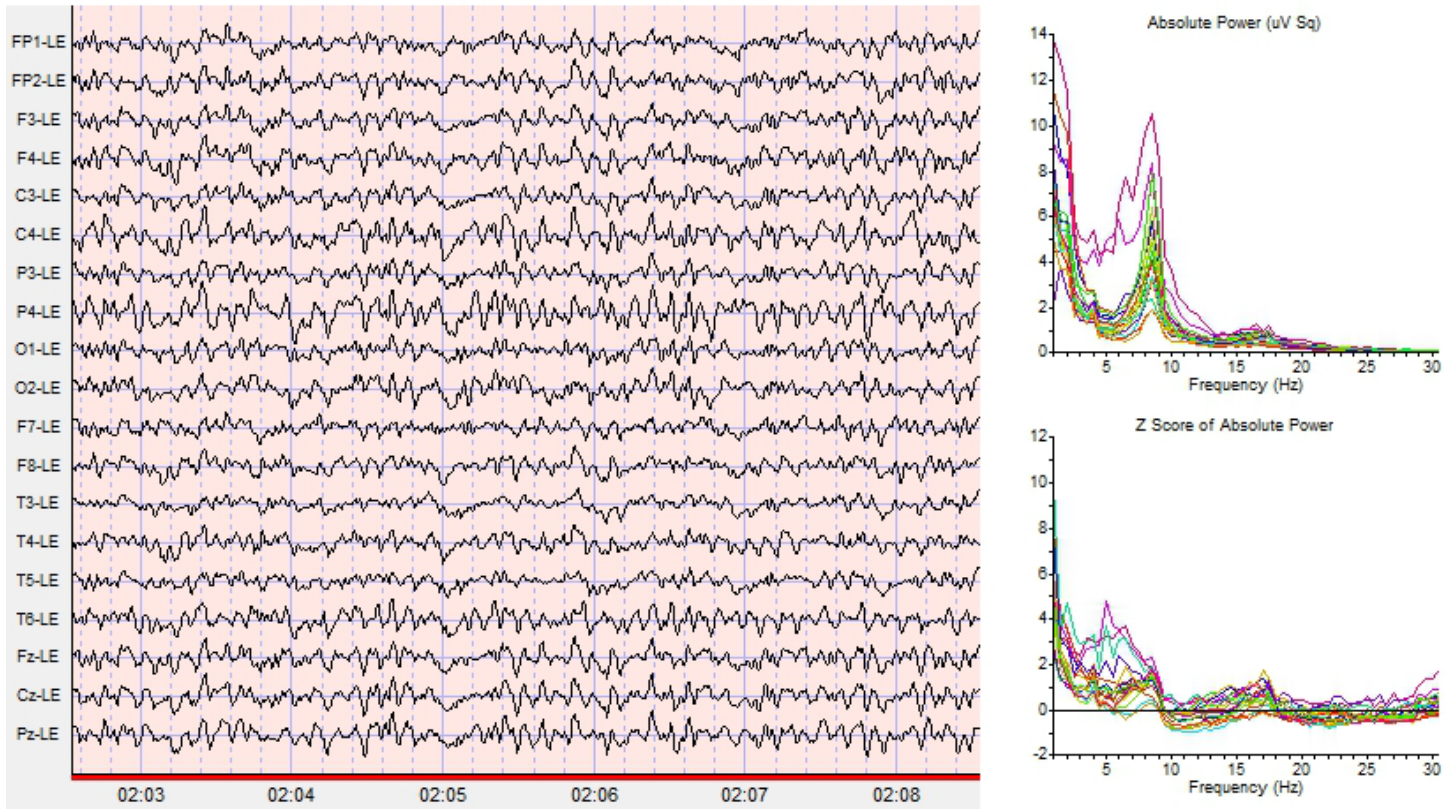


Example of Laplacian EEG and Absolute Power - Eyes Closed Condition

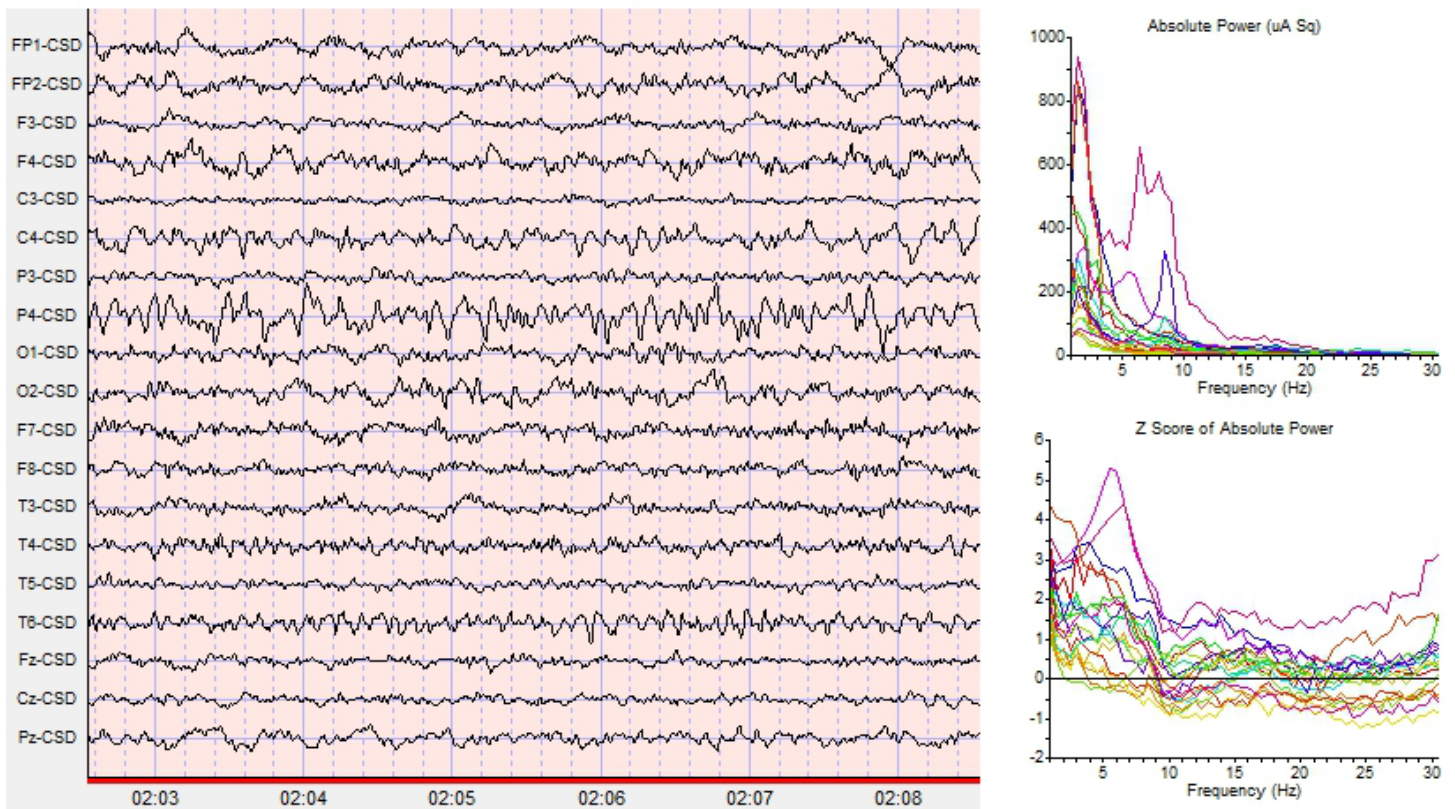

Electrical Neuroimaging


Linking a patient's symptoms and complaints to functional systems in the brain is important in evaluating the health and efficiency of cognitive and perceptual functions. The electrical rhythms in the EEG arise from many sources but approximately 50% of the power arises directly beneath each recording electrode. Electrical Neuroimaging uses a mathematical method called an "Inverse Solution" to accurately estimate the sources of the scalp EEG (Pascual-Marqui et al, 1994; Pascual-Marqui, 1999). Below is a Brodmann map of anatomical brain regions that lie near to each 10/20 scalp electrode with associated functions as evidenced by fMRI, EEG/MEG and PET Neuroimaging methods.

Symptoms, Electrodes & Brodmann Areas

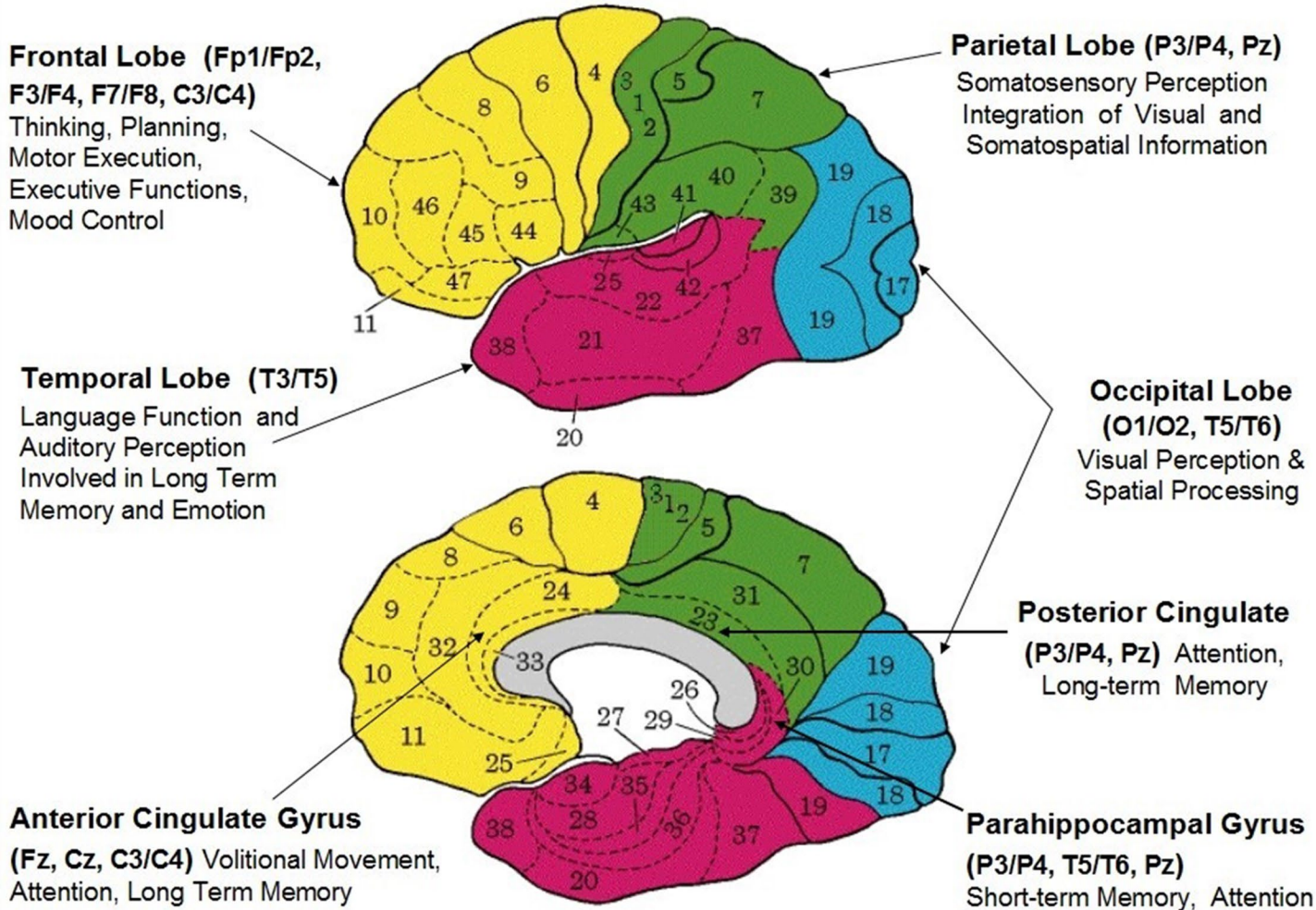


Significant Brain Lobe Z-Scores

FRONTAL BRODMANN AREAS

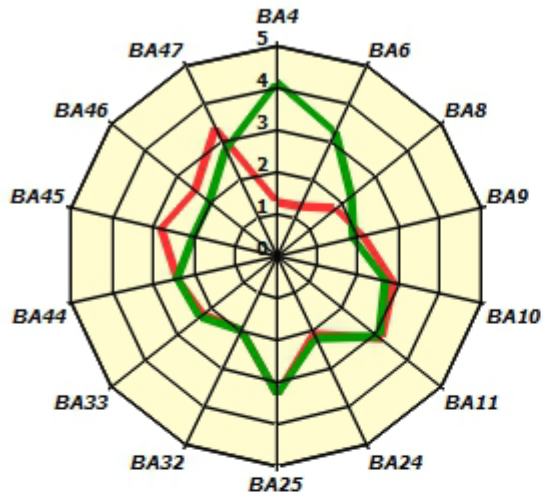


TEMPORAL BRODMANN AREAS

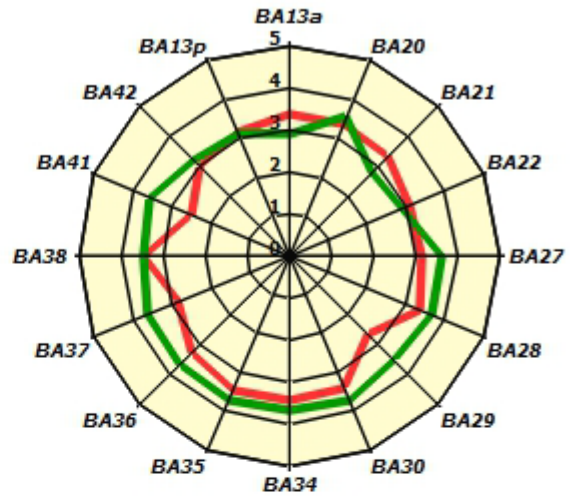


■ LEFT ■ RIGHT

PARIETAL BRODMANN AREAS

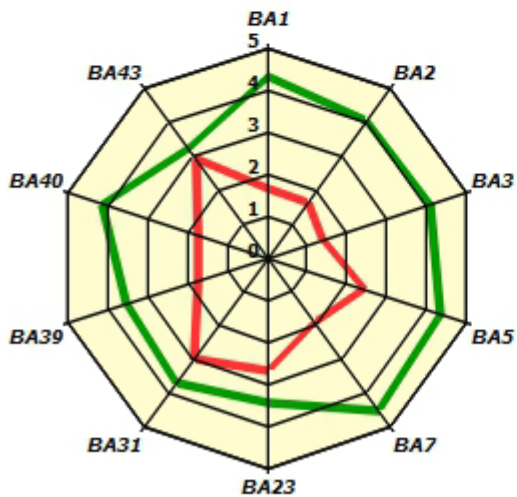

OCCIPITAL BRODMANN AREAS

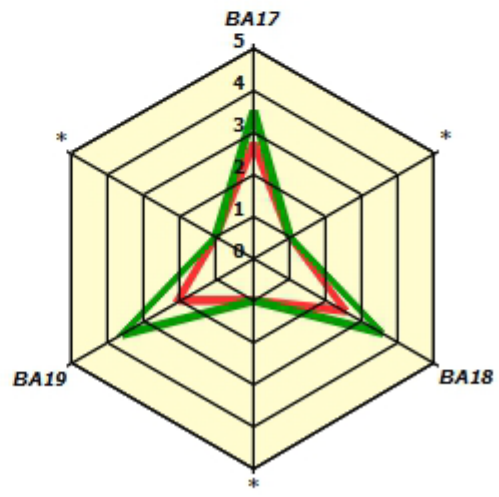


Significant LORETA NeuroImaging at 1 Hz - Brodmann Area 30

Value = 4.10301 - 1 Hz




© 2016-2023 Applied Neuroscience, Inc.

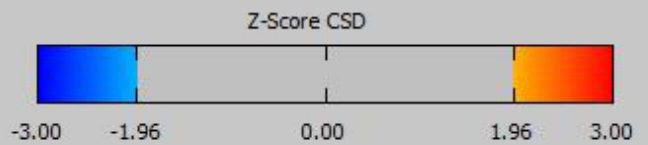


Significant LORETA NeuroImaging at 2 Hz - Brodmann Area 7

Value = 3.84424 - 2 Hz

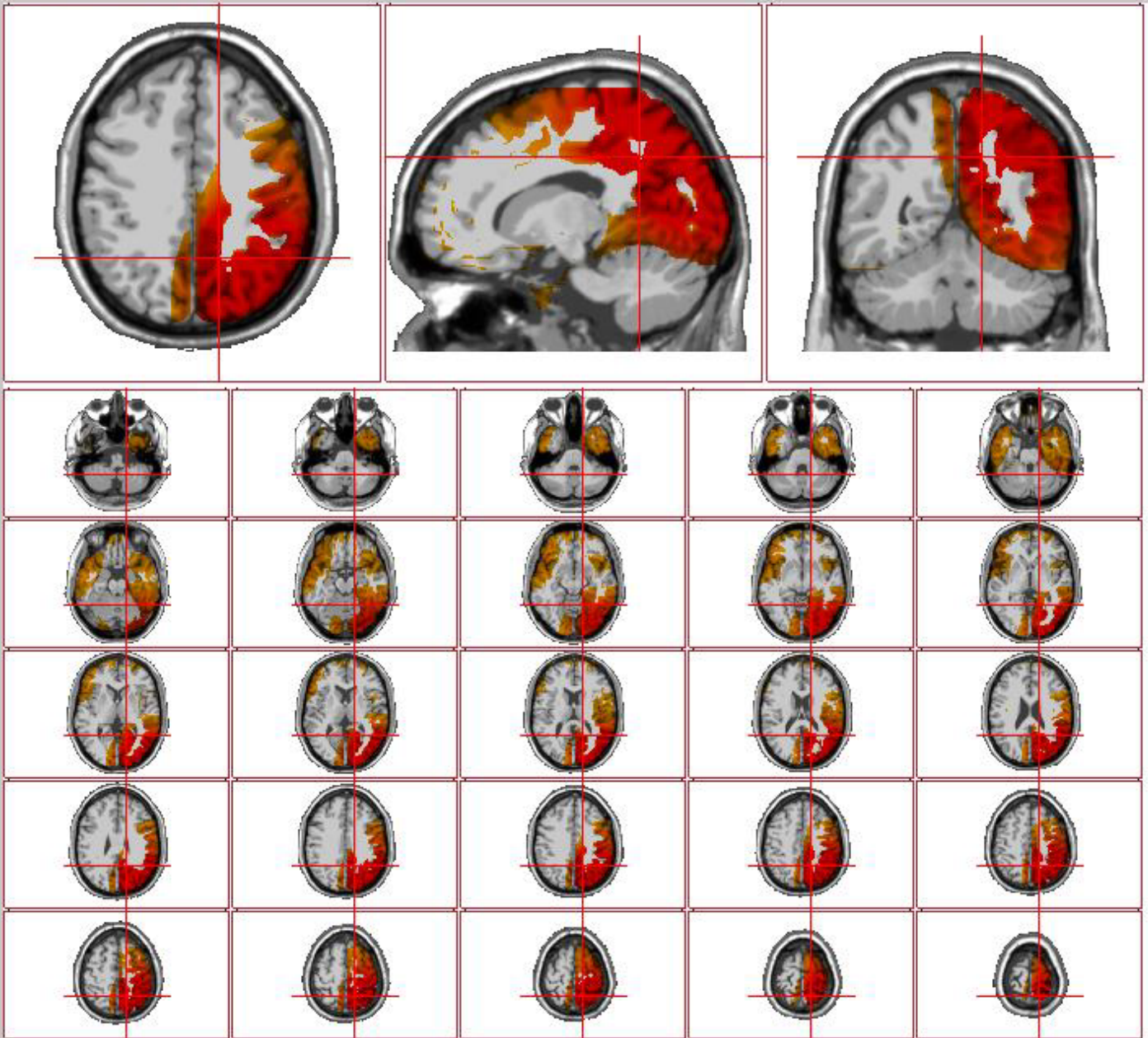


© 2016-2023 Applied Neuroscience, Inc.

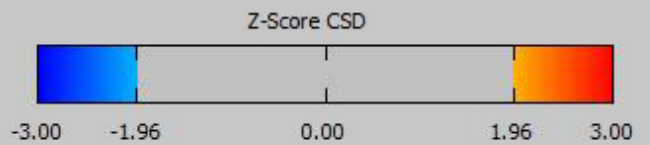


Significant LORETA NeuroImaging at 4 Hz - Brodmann Area 7

Value = 4.38533 - 4 Hz

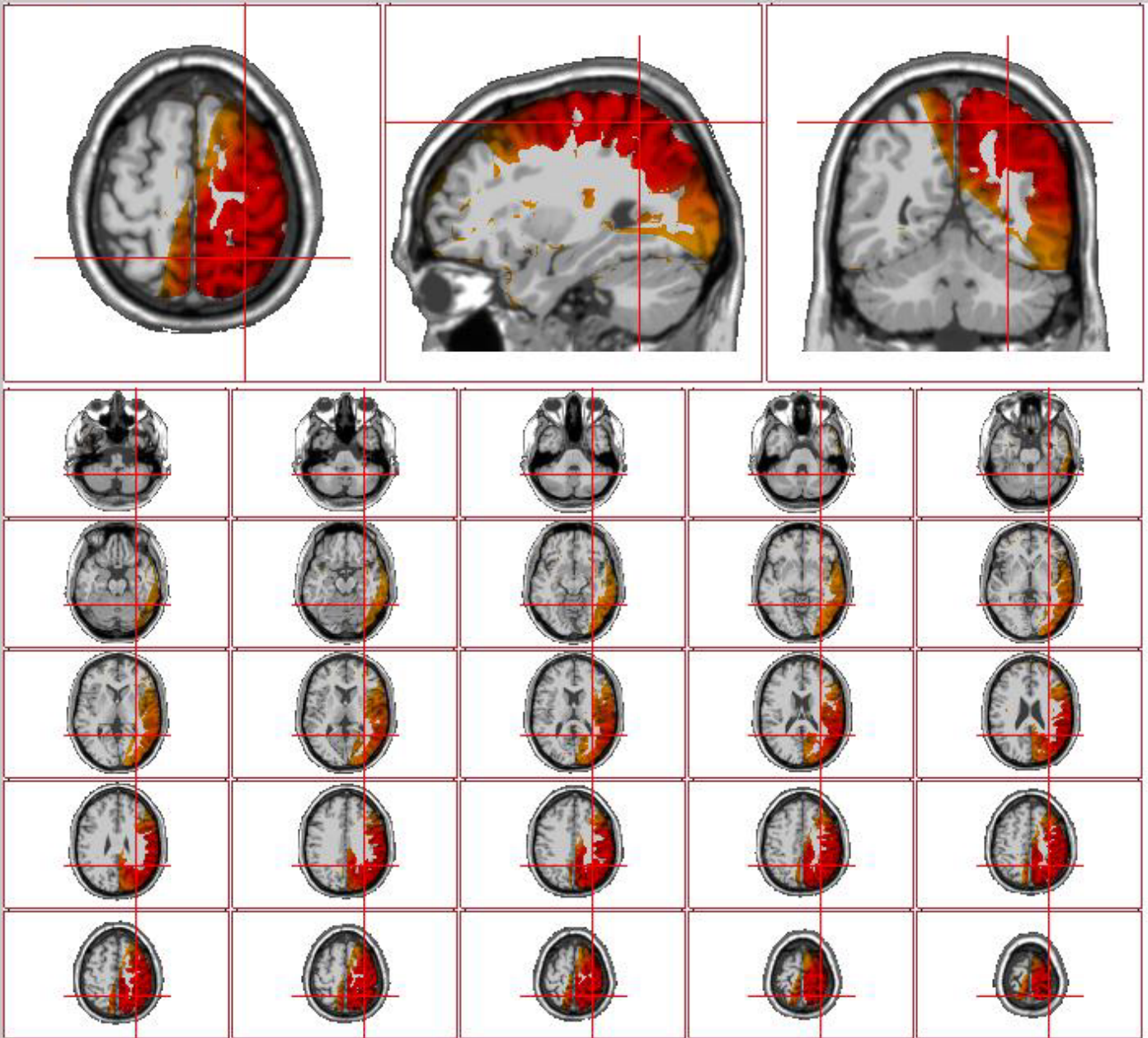


© 2016-2023 Applied Neuroscience, Inc.

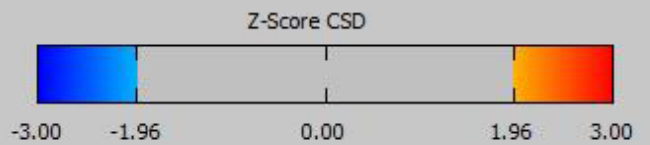

Significant LORETA NeuroImaging at 5 Hz - Brodmann Area 7


Value = 4.83668 - 5 Hz

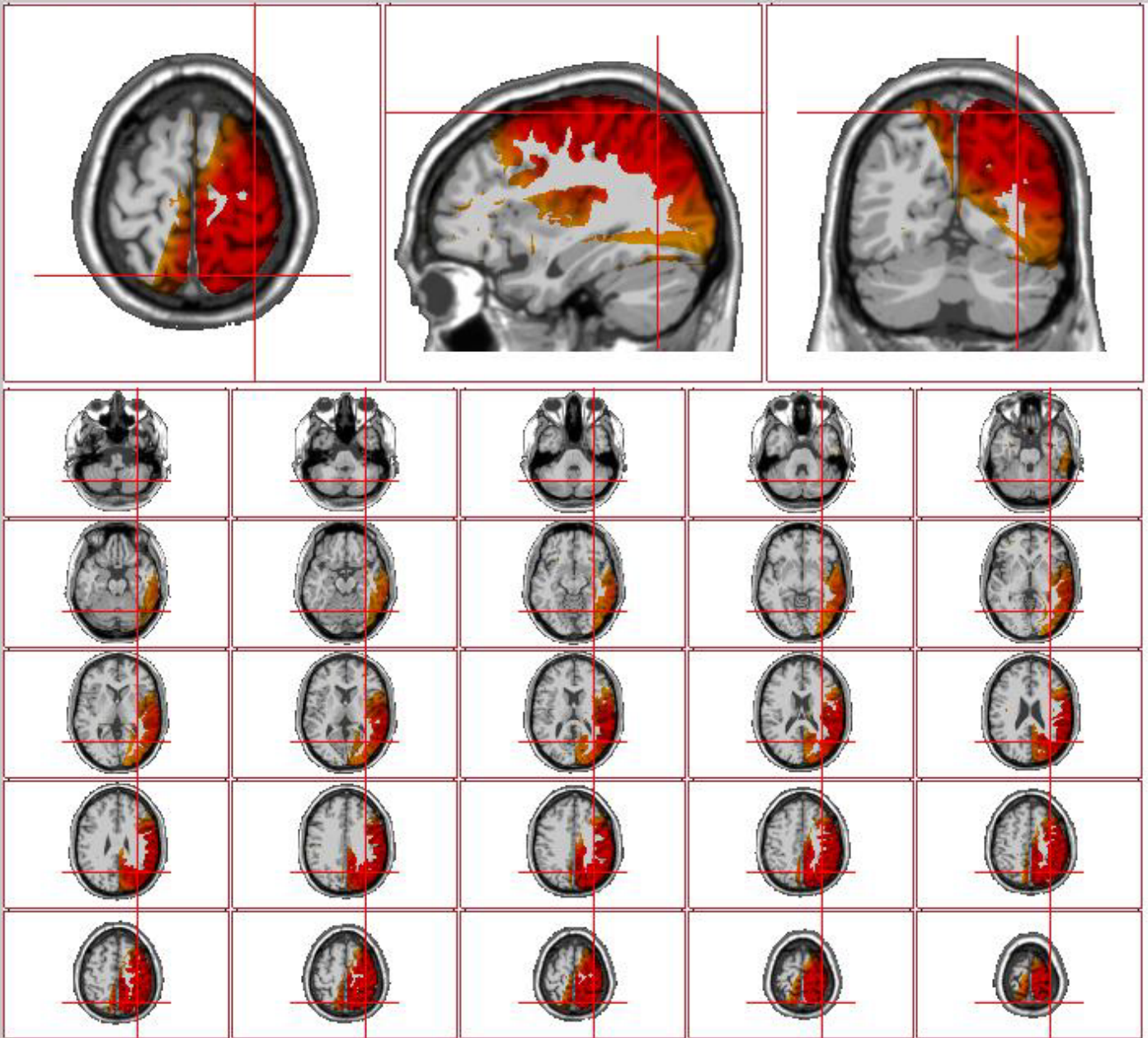


© 2016-2023 Applied Neuroscience, Inc.


Significant LORETA NeuroImaging at 6 Hz - Brodmann Area 1


Value = 4.50427 - 6 Hz

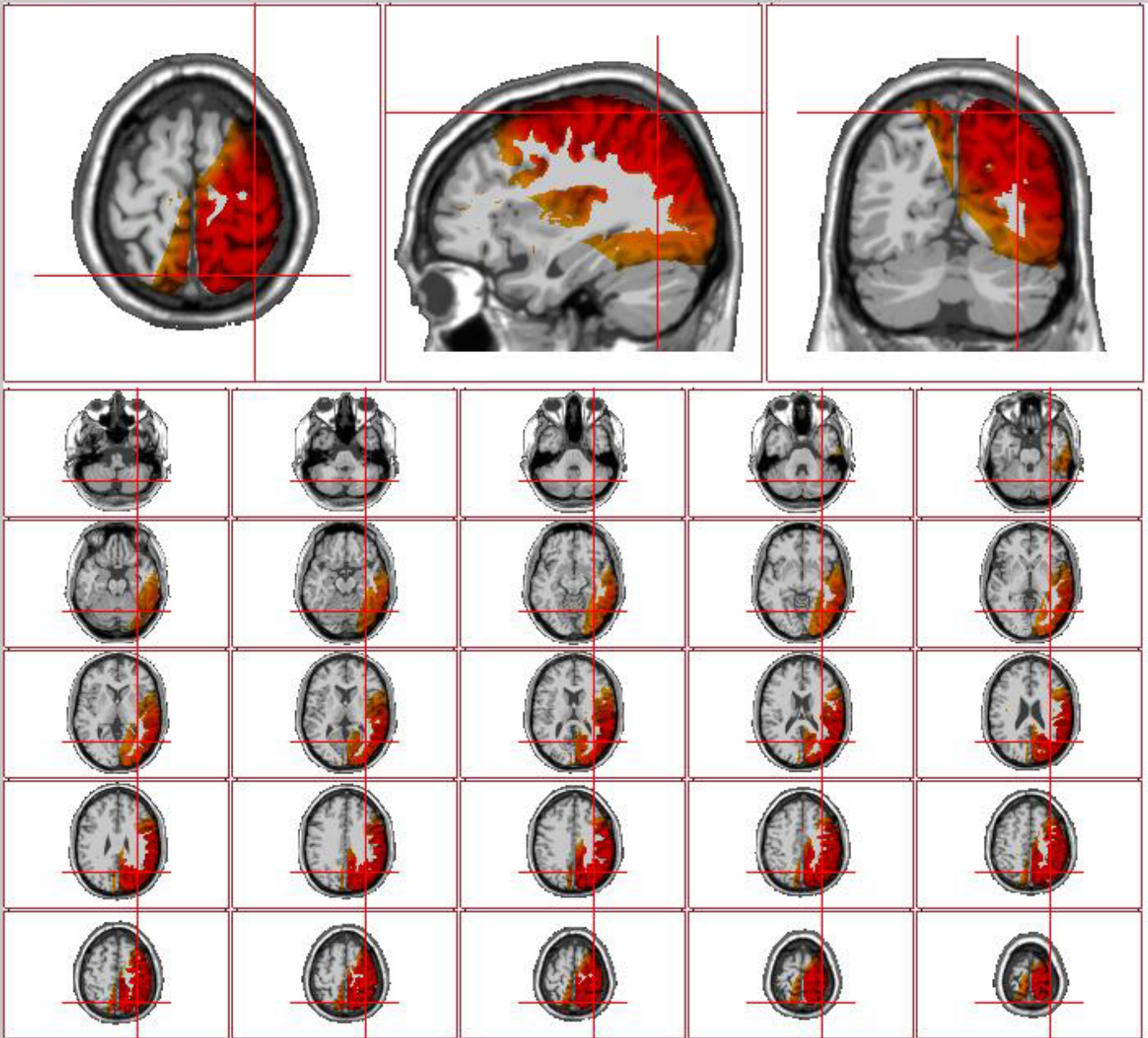



© 2016-2023 Applied Neuroscience, Inc.

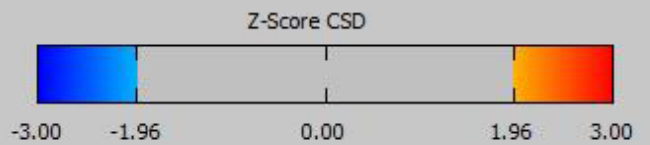

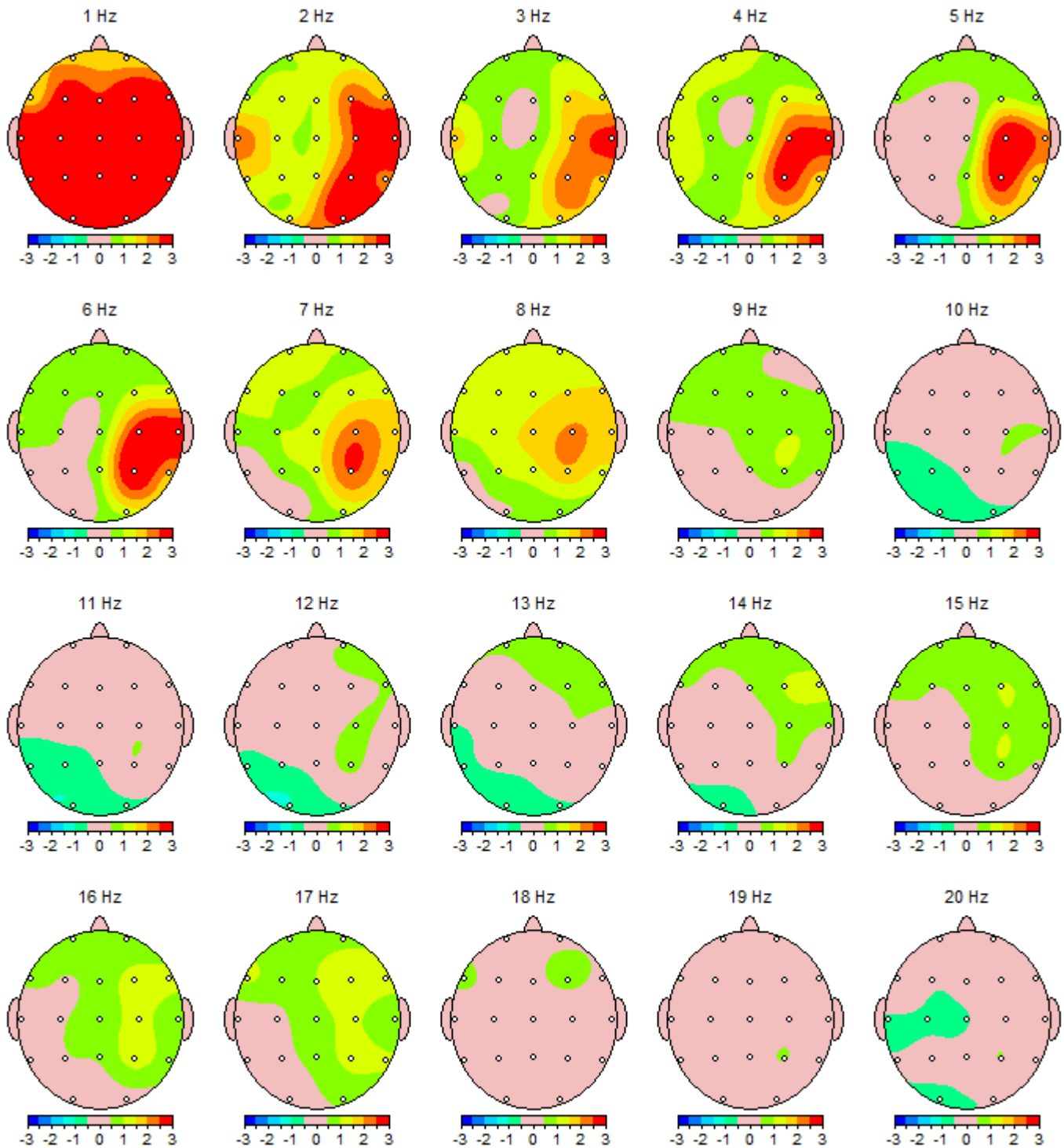
Linked-Ears Absolute Power (1-30 Hz) QEEG Maps

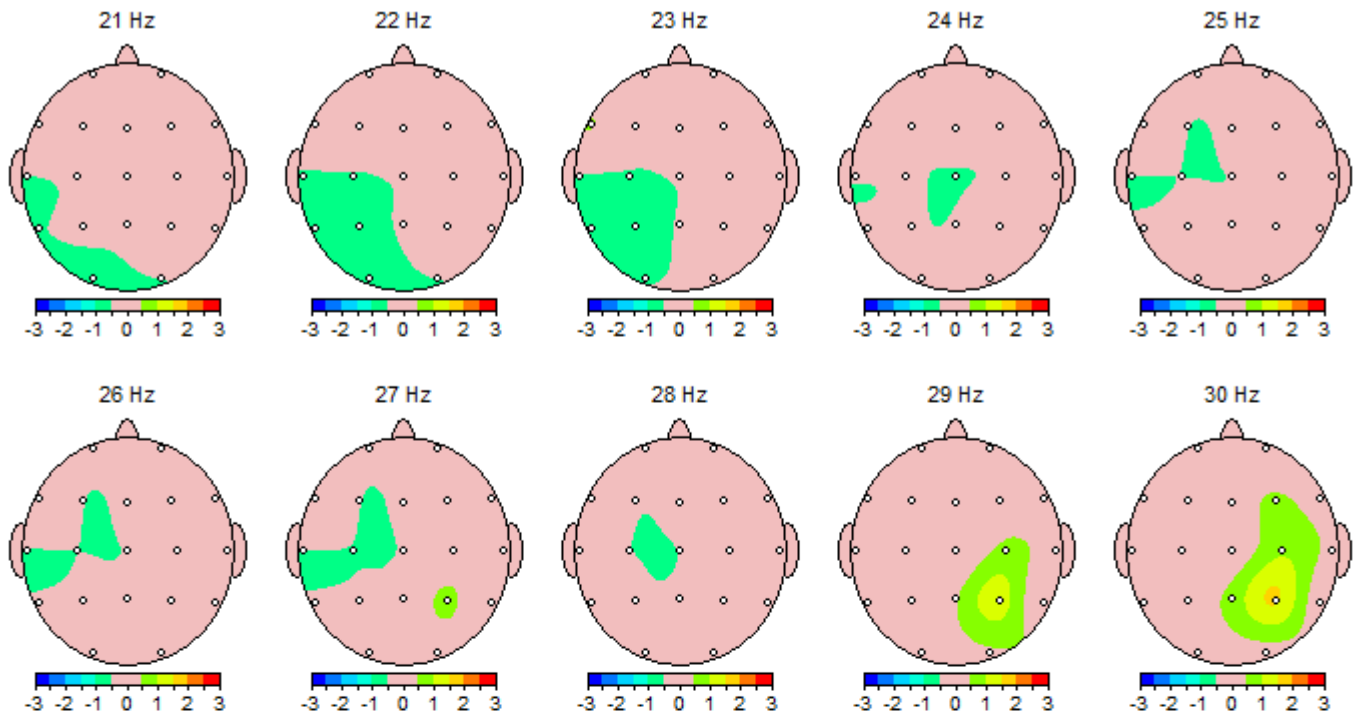Linked-Ears Absolute Power (1-30 Hz) QEEG Maps

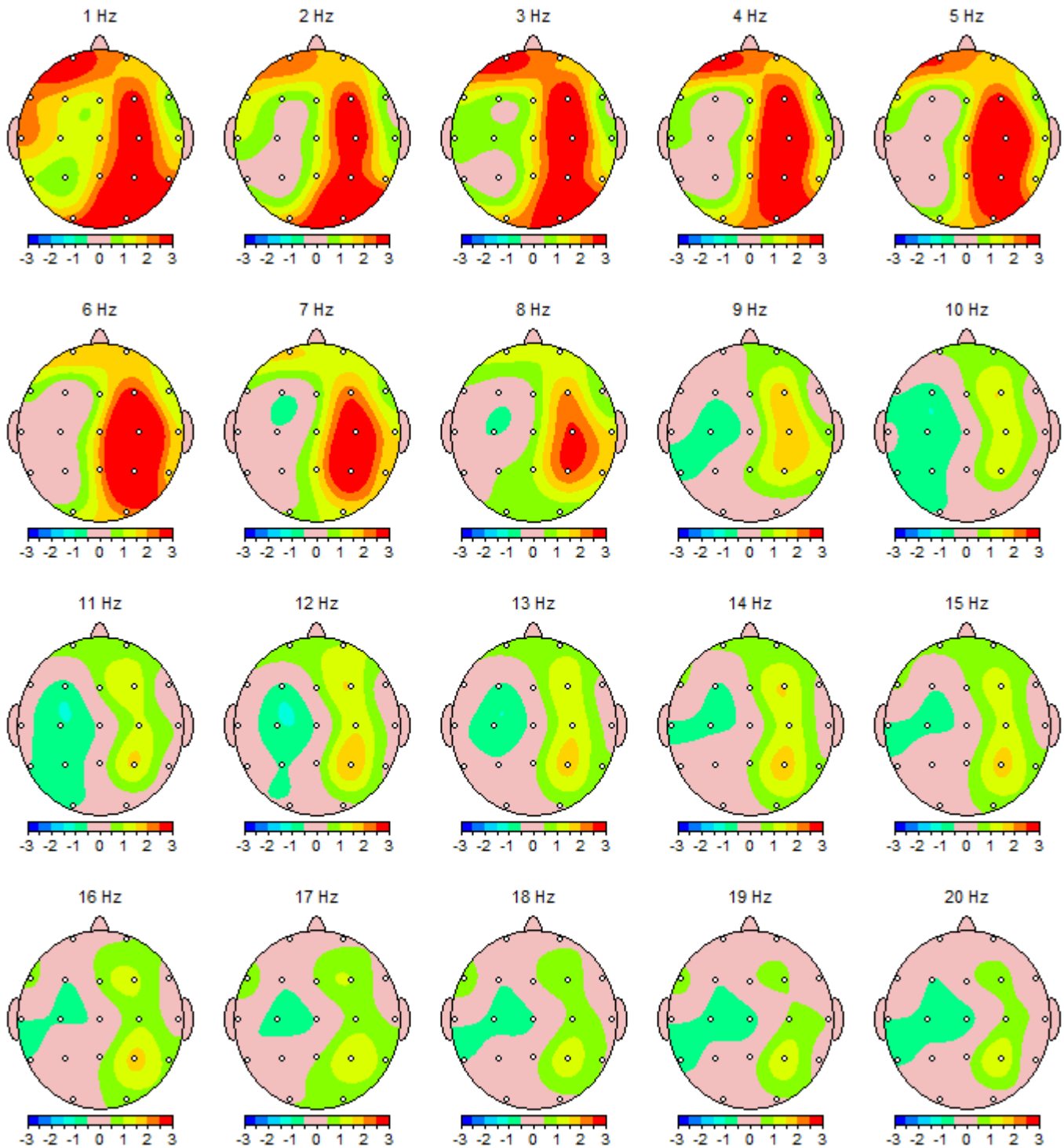
Laplacian Absolute Power (1-30 Hz) QEEG Maps

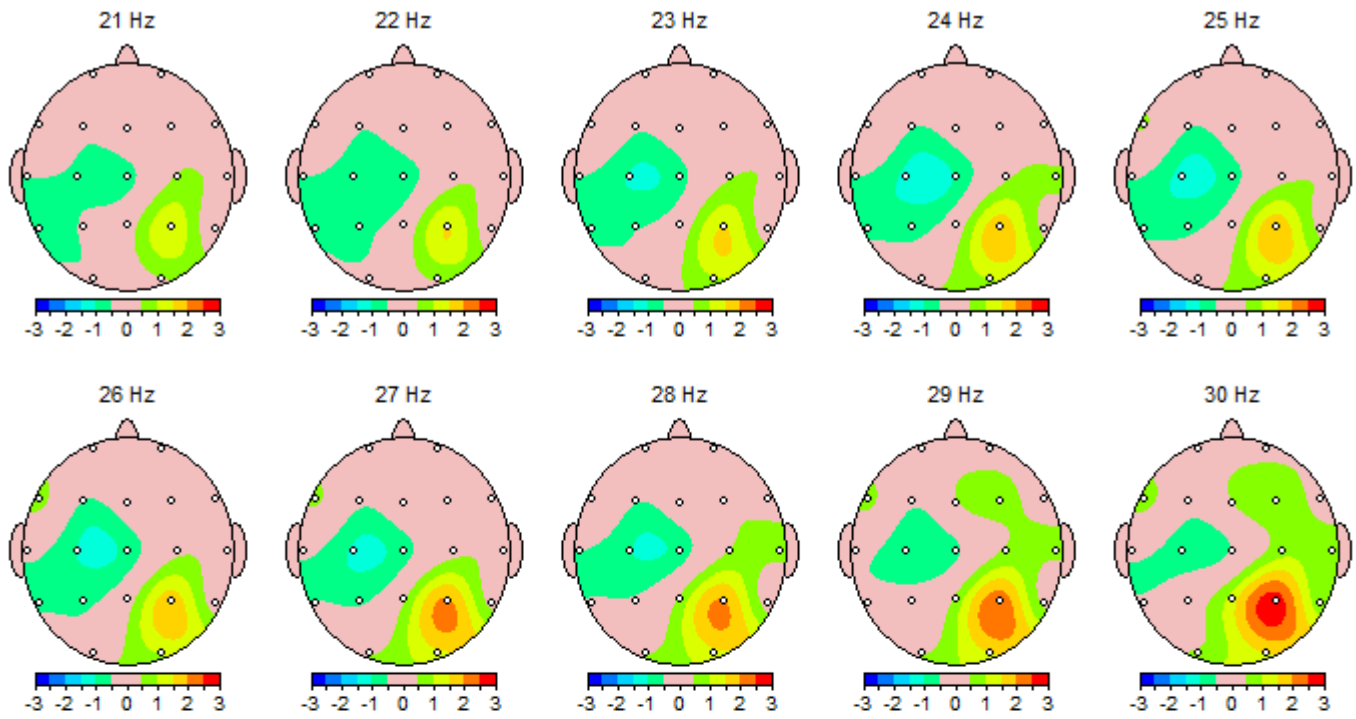Laplacian Absolute Power (1-30 Hz) QEEG Maps


Significant Z-Scores Amplitude Asymmetry Connection Maps

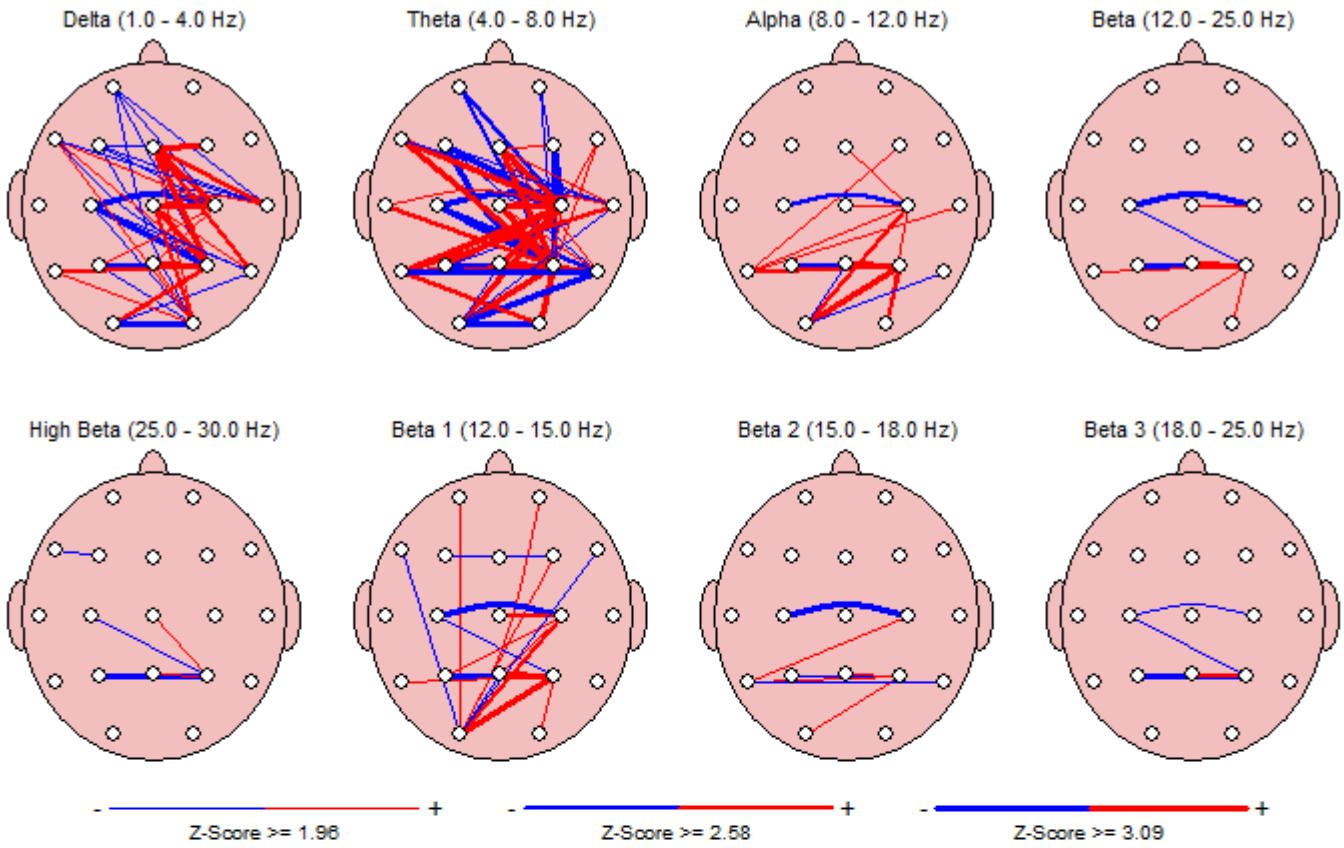

Significant Z-Scores Coherence Connection Maps

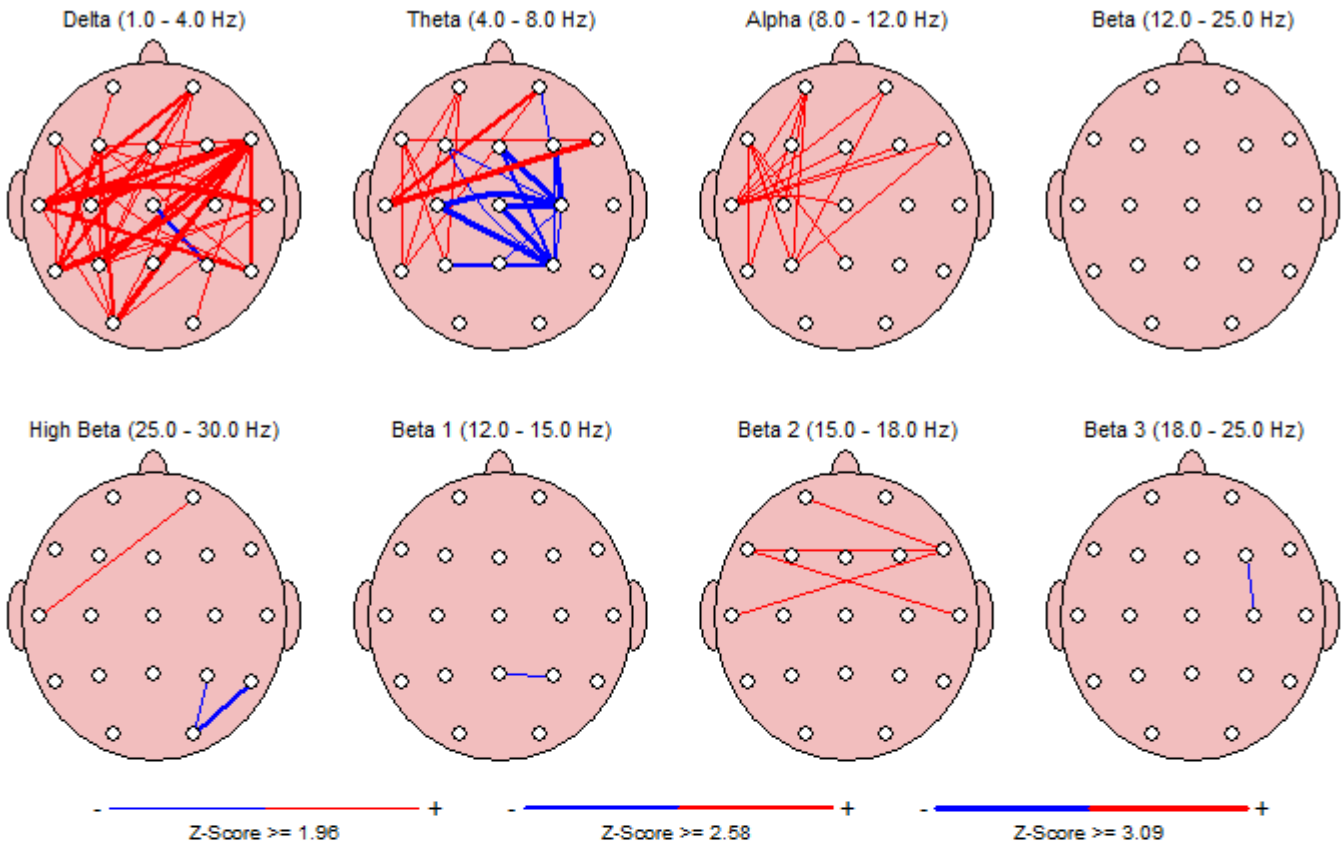

Significant Z-Scores Phase Lag Connection Maps

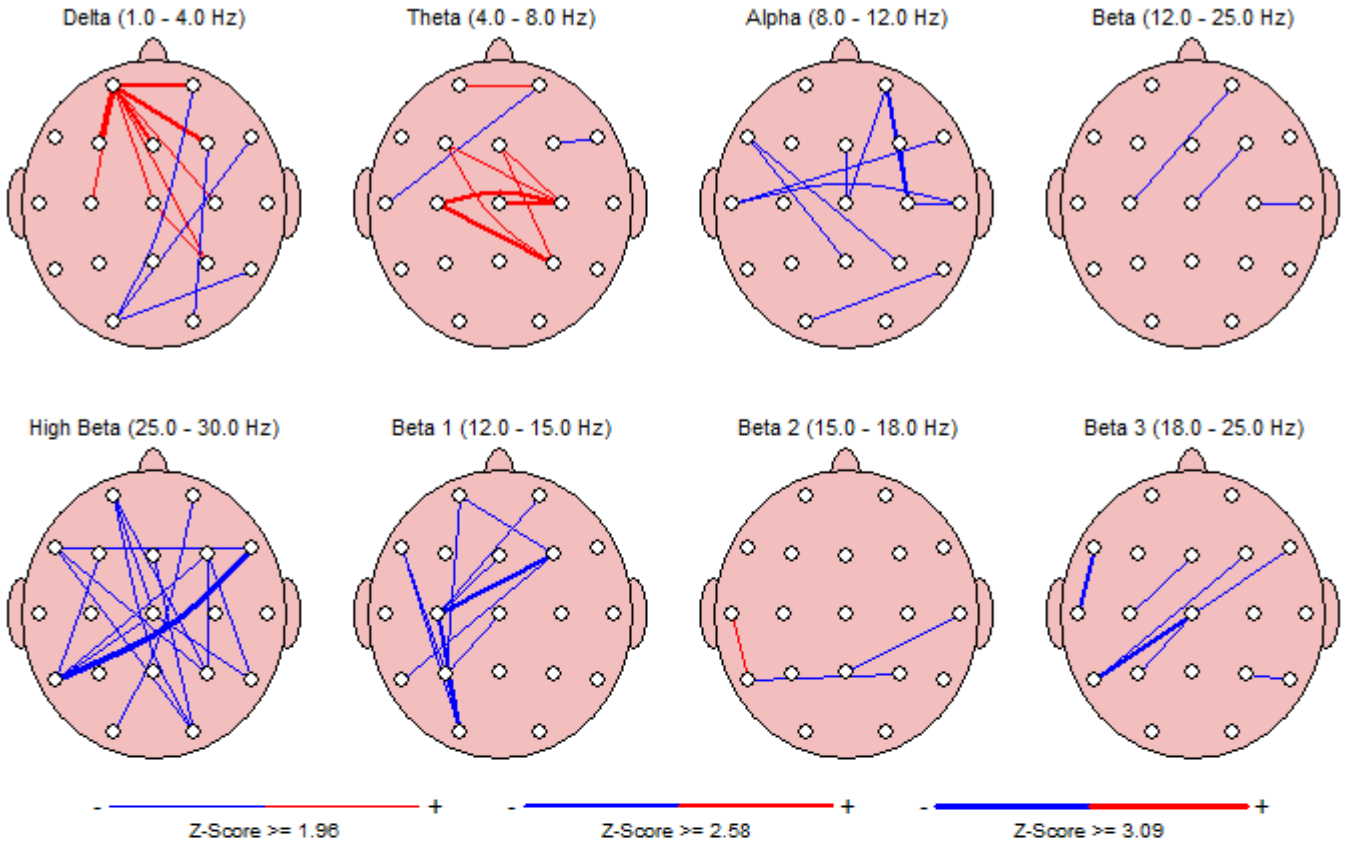


Traumatic Brain Injury Discriminant Analysis

TBI DISCRIMINANT SCORE = 1.15

TBI PROBABILITY INDEX = 99.5%


The TBI Probability Index is the subject's probability of membership in the mild traumatic brain injury population. (see Thatcher et al, EEG and Clin. Neurophysiol., 73: 93-106, 1989.)

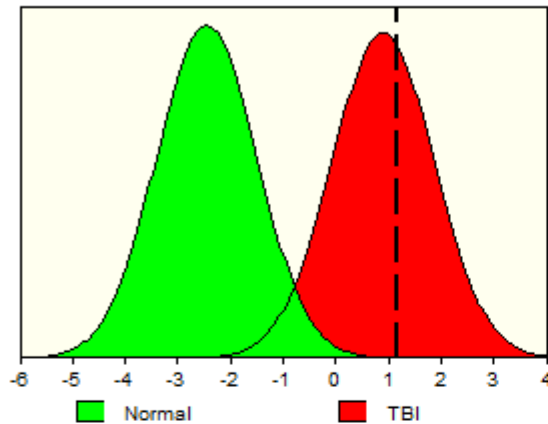

			RAW	Z
FP1-F3	COH	Theta	85.64	0.41
T3-T5	COH	Beta	70.27	1.53
C3-P3	COH	Beta	87.89	1.77
FP2-F4	PHA	Beta	0.26	-0.33
F3-F4	PHA	Beta	0.40	0.18
F4-T6	AMP	Alpha	-7.22	0.25
F8-T6	AMP	Alpha	-60.32	0.09
F4-T6	AMP	Beta	30.19	0.73
F8-T6	AMP	Beta	-20.24	0.48
F3-O1	AMP	Alpha	34.05	1.82
F4-O2	AMP	Alpha	14.84	1.36
F7-O1	AMP	Alpha	6.31	-2.33
F4-O2	AMP	Beta	32.07	1.49
P3	RP	Alpha	38.19	-1.05
P4	RP	Alpha	38.15	-1.12
O1	RP	Alpha	34.11	-1.56
O2	RP	Alpha	27.93	-2.09
T4	RP	Alpha	30.86	-1.08
T5	RP	Alpha	30.03	-1.62
T6	RP	Alpha	39.66	-1.07

TBI SEVERITY INDEX = 3.60

This severity score places the patient in the MODERATE range of severity.

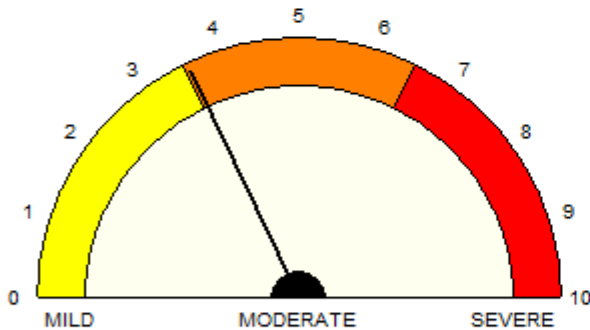

			RAW	Z
FP1-C3	COH	Delta	58.70	0.58
FP1-FP2	COH	Theta	72.74	-1.94
O1-F7	COH	Alpha	36.72	0.74
O2-T6	COH	Alpha	79.58	-0.12
P3-O1	COH	Beta	81.97	0.79
FP1-T3	PHA	Theta	4.88	0.76
T3-T4	PHA	Theta	-7.38	-0.24
O1-F7	PHA	Alpha	-7.38	-0.29
F7-F8	PHA	Alpha	0.89	-0.68
T5-T6	PHA	Beta	2.02	0.00
C3-F7	AMP	Delta	29.06	-0.59
FP2-F4	AMP	Delta	-29.04	-0.57
C4-F8	AMP	Delta	58.80	0.63
O1-O2	AMP	Theta	-72.90	-4.90
P3-F7	AMP	Alpha	33.91	-1.78
FP2-P4	AMP	Alpha	-100.31	-0.55

The TBI Severity Index is an estimate of the neurological severity of injury. (see Thatcher et al, J. Neuropsychiatry and Clinical Neuroscience, 13(1): 77-87, 2001.)

***Statement of Indications of Use:**

The Discriminant Analysis and Severity Index are to be used by experienced and qualified professionals for the post-hoc statistical evaluation of the human electroencephalogram (EEG). The Discriminant Analysis and Severity Index are to be viewed as an adjunct to the evaluation of the patient, and they do not serve as a primary basis for a diagnosis. Warning: Inclusion criteria of a history of traumatic brain injury and greater than 13 years of age must be adhered to.

Quality Control of EEG Collection

Record Length: 03:48

Edit Length: 02:19

Reliability:

	Split Half	Test Retest
Average	0.99	0.98
FP1	0.99	1.00
FP2	0.99	0.94
F3	1.00	0.99
F4	0.97	0.99
C3	1.00	0.98
C4	0.99	0.99
P3	1.00	0.97
P4	1.00	0.98
O1	0.98	0.99
O2	0.97	0.98
F7	0.99	0.96
F8	0.98	0.99
T3	1.00	0.98
T4	0.98	0.99
T5	0.99	0.97
T6	0.98	0.96
Fz	1.00	0.99
Cz	1.00	0.97
Pz	1.00	0.94

Sampling Rate: 128

Collection Hardware: Lexicor (High Pass Off)

An Addendum to NeuroGuide QEEG Report

Important Disclaimer: QEEG tests are ancillary tests that are not intended to provide a diagnosis by themselves, but are used to evaluate the nature and severity of deregulation in the brain such as in mild traumatic brain injury (MTBI). The QEEG tests provide a quantitative assessment of areas of brain dysfunction and information on impaired conduction and connectivity between different regional neural networks in the brain. The assessment of impaired connectivity is based on abnormal measurements of Coherence and Phase. The TBI Discriminant and Concussion Index do not provide a diagnosis for MTBI but only information on the presence of a pattern in the EEG that is often found in patients with a history of mild traumatic brain injury. The TBI Discriminant and Concussion Index also provide information about connectivity and excitability of brain regions. The TBI Discriminant and Concussion Index are to be used only on patients with a clinical history and symptoms of a Traumatic Brain Injury and Post Concussion syndrome. The diagnosis of MTBI is a clinical one and is not based on any one test. A diagnosis is performed by the clinician, who integrates the medical history, clinical symptoms, neurocognitive tests with the abovementioned brain function tests as well as other information to render a diagnosis. The information on impaired brain connectivity is derived primarily from abnormal measurements of Coherence and Phase. Assessments of regional abnormality rely also on abnormal amplitude (power) distribution across the spectrum of EEG frequencies as compared to the normative database.

Artifact Rejection: NeuroGuide uses the standard deletion of artifact method to only select artifact free EEG data for analyses. The entire EEG record must be viewed by clicking end and page down and page up and home and by arrow keys and by moving the wiper at the bottom of the screen. A careful visual examination of the EEG record is necessary to detect epilepsy and gross pathology as well as to identify artifacts. The goal is to avoid selecting any artifact and instead to only select artifact free segments of EEG. There are three methods of obtaining Artifact Free Selections: 1- Manual Selections are obtained by pressing the left mouse button and dragging to select, press right mouse button and drag to erase; 2- Artifact Free Template Matching; and 3- Z Score Artifact Free Selections. All three methods can be used and manual selection takes priority over all methods of artifact free selection. That is, left and right mouse button dragging will override all other methods. View the Length of EEG Selections in seconds and View the dynamic Reliability Measures of the EEG Selections. For Manual Selections of Artifact Free EEG Depress the left mouse button and drag it over the sections of EEG that do not contain eye movement or muscle or drowsiness or head movement or any other type of artifact. Select at least 60 seconds of artifact free EEG data as shown in the Edit Time counter (upper left of screen). If a mistake is made, then right mouse click and drag over the EEG traces to erase a selection. View the Test Re-Test reliability which must be at least 0.90. Scan the EEG record and select real and valid EEG and avoid selecting artifact. Splice discontinuities are removed by filtering and exercises to prove no distortion due to splicing are available in the Handbook of QEEG and EEG Biofeedback. Pattern recognition routines are used to identify likely eye movement (EOG), drowsiness and muscle (EMG) artifact in the record and thereby mark these suspected segments and disallow them to be included in subsequent analyses. The pattern recognition routines are based on physics and physiology of artifact. For example, all electrical sources decrement with distance and in the case of eye movement detection is by the presence of an electrical field gradient in the delta frequency band from $Fp1/2 > F3/4 > C3/4$ and/or 120 degrees or higher of inverse phase between F7 and F8. EMG electrical gradients at > 10 Hz from $T3/4 > C3/4$ and/or $Fp1/2 > F3/4 > C3/4$ and/or $O1/2 > P3/4$. Drowsiness occurs when the locus coeruleus reduces inhibition on the hypothalamic sleep centers resulting in 2 - 4 Hz action potential bursting that projects to the ventral posterior thalamic relay nuclei. Drowsiness pattern detection involves elevated

slow waves in the EEG maximal in Cz and Fz as well as alpha slowing. NeuroGuide does not use any regression methods to allegedly remove artifact such as ICA/PCA or Blind Source or unpublished methods like SARA that distort Phase and Coherence and other aspects of the Power Spectrum. Details and tutorials demonstrating how the ICA and regression methods distort Phase and Coherence are available at:

https://www.appliedneuroscience.com/PDFs/Tutorial_Adulteration_Phase_Relations_when_using_ICA.pdf

Split Half and Test Re-Test Reliability: Split-Half (SH) reliability is the ratio of variance between the even and odd seconds of the time series of selected digital EEG (variance = sum of the square of the deviation of each time point from the mean of the time points). Examine the average reliability and the reliability of each channel as you increase the length of the sample and manually select different segments. Selection of artifact free EEG should have a reliability > 0.95 and a sample length of edited EEG > 60 seconds. Test Re-Test (TRT) reliability is the ratio of variance between the first half vs. the second half of the selected EEG segments (variance = sum of the square of the deviation of each time point from the mean of the time points). Test Re-Test reliability > 0.90 and a sample length of edited EEG > 60 seconds is commonly published in the scientific literature. Test Re-Test reliability is an excellent statistic to compare Brain state changes such as drowsiness as well as the consistency of a measure independent of changes in brain state.

NeuroGuide Normative Database Description: The NeuroGuide normative database in versions 1.0 to 2.4.6 included a total of 678 carefully screened individual subjects ranging in age from 2 months to 82 years. NG 2.6.8 involved the addition of 49 adult subjects ranging in age from 18.3 years to 72.6 years resulting in a normative database of 727 subjects. The inclusion/exclusion criteria, demographics, neuropsychological tests, Gaussian distribution tests and cross-validation tests are described in several peer reviewed publications (Thatcher et al, 1983; 1987; 2003). Two year means were computed using a sliding average with 6 month overlap of subjects. This produced a stable and higher age resolution normative database with a total of 21 different age groups. The 21 age groups and age ranges and number of subjects per age group is shown in the bar graph in Appendix F figure 2 in the NeuroGuide Manual (click Help > NeuroGuide Help).

www.appliedneuroscience.com

and clicking on the webpage Articles & Links > Articles > Article #5.

Normative Database for Different Montages: The raw digital data from the same group of normal subjects is analyzed using different montages such as Average Reference, Laplacian current source density, a common reference based on all 19 channels of the 10/20 system and standard clinical bipolar montages (e.g., longitudinal, circular, transverse). Users can create any montage that they wish and there will be a normative reference database comparison available for both eyes closed and eyes open conditions.

Age Range of LORETA Normative Databases: The LORETA current density and source correlation norms use the same subjects as are used for the surface EEG norms and the age range is 2 months to 82 years. The computational details of the LORETA current density norms are published at: Thatcher, R.W., North, D., Biver, C. EEG inverse solutions and parametric vs. non-parametric statistics of Low Resolution Electromagnetic Tomography (LORETA). Clin. EEG and Neuroscience, 36(1): 1-9, 2005 and Thatcher, R.W., North, D., Biver, C. Evaluation and Validity of a LORETA normative EEG database. Clin. EEG and Neuroscience, 2005, 36(2): 116-122. Copies of these publications are available to download from www.appliedneuroscience.com

by clicking Articles & Links > Articles > Numbers 11 and 12. The computational details of the LORETA source correlation norms are in the NeuroGuide Manual, click Help > NeuroGuide Help > Appendix-G.

Implementation of LORETA Measurement in NeuroGuide: The Key Institute's LORETA equations and the LORETA viewer (Pascual-Marqui et al, 1994; Pascual-Marqui, 1999) can be launched by a single mouse click in the NeuroGuide window. NeuroGuide exports frequency domain and time domain edits of 19 channel x 256 point digital EEG in microvolts (or μV^2) in the Lexicor electrode order as the standard input to the Key Institute T-Matrix. Rows are 256 microvolt time points and the columns are 19 channels at a sample rate of 128 thus producing 0.5 Hz resolution from 1 to 30 Hz. 1 Hz increments in the LORETA viewer are computed as the sum of adjacent 0.5 Hz bins and thus the 'Time Frame' control in the LORETA Viewer is frequency from 1 to 30 Hz. See Pascual-Marqui RD, Michel CM, Lehmann D., 1994. Low resolution electromagnetic tomography: a new method for localizing electrical activity in the brain. *International J. of Psychophysiology*, 18:49-65. For computational details see: Pascual-Marqui. R.D., 1999. Review of Methods for Solving the EEG Inverse Problem. *International J. of Bioelectromagnetism*, 1(1): 75-86. Pascual-Marqui, R.D., 2004. The Key Institute's free software and documentation was downloaded from www.unizh.ch/keyinst/NewLORETA/Software/Software.htm

Amplifier Matching Importance: This stems from the fact that amplifiers have different frequency gain characteristics. The matching of amplifiers to the NeuroGuide database amplifier was done by injecting microvolt calibration signals of different amplitudes and frequencies into the input of the respective EEG machines and then computing correction curves to exactly match the amplifier characteristics of the norms and discriminant functions. The units of comparison are in microvolts and a match within 3% is generally achieved. The NeuroGuide research team double checked the amplifier match by computing FFT and digital spectral analyses on calibration signals used to acquire the norms with the calibration signals used to evaluate a given manufacturers amplifiers.

Historical QEEG Normative Databases Scientific Standards: A review of the history of QEEG normative databases was published in Thatcher, R.W. and Lubar, J.F. History of the scientific standards of QEEG normative databases. In: *Introduction to QEEG and Neurofeedback: Advanced Theory and Applications*, T. Budzinsky, H. Budzinsky, J. Evans and A. Abarbanel (eds)., Academic Press, San Diego, CA, 2008. A copy of the publication can be downloaded at:

https://www.appliedneuroscience.com/PDFs/History_of_QEEG_Databases.pdf

QEEG Normative Database Publications and Validations

Bosch-Bayard J, Valdes-Sosa P, Virues-Alba T, Aubert-Vazquez E, John ER, Harmony T, Riera-Diaz J, Trujillo-Barreto N. (2001). 3D statistical parametric mapping of EEG source spectra by means of variable resolution electromagnetic tomography (VARETA). *Clin Electroencephalogr.*, 32(2):47-61.